

STEADFAST RESEARCH GROUP

Label & Packaging Guidelines

Sizing, layout, materials, and print specification for peptide vial labels and unit packaging. Version 1.0 · July 2026.

FOR RESEARCH USE ONLY · NOT FOR HUMAN CONSUMPTION

Know the glass before the label.

SRG products ship in USP Type I borosilicate glass vials with 13 mm crimp necks, flip-off caps in olive, and lyophilized contents. Label size is driven by two numbers: the **body circumference** (label width) and the **straight sidewall height** between shoulder and heel (label height).


VIAL	BODY (Ø × H)	CIRCUMFERENCE	STRAIGHT WALL	RECOMMENDED LABEL
2 mL	Ø16 × 35 mm	50.3 mm	≈22 mm	47 × 19 mm
3 mL (standard)	Ø22 × 47 mm	69.1 mm	≈30 mm	63.5 × 25.4 mm (2.5" × 1")
5 mL	Ø24 × 52 mm	75.4 mm	≈34 mm	70 × 30 mm
10 mL (future)	Ø25 × 54 mm	78.5 mm	≈36 mm	73 × 32 mm

Wrap rules

- **Coverage:** label width ≈90% of circumference. Leave a 5–8 mm vertical gap of bare glass - do not overlap label onto itself; overlap lifts at freezer temperatures.
- **Vertical position:** center the label on the straight wall with ≥2 mm clearance from shoulder and heel. Label applied over the curve will flag (peel at edges).
- **Application:** machine-applied wrap, seam at the back, artwork centered opposite the seam.

Confirm before ordering: vial dimensions vary ±1 mm between glass suppliers. Measure a physical sample of the actual vial lot before sending label artwork to print.

The 3 mL label, dimensioned.

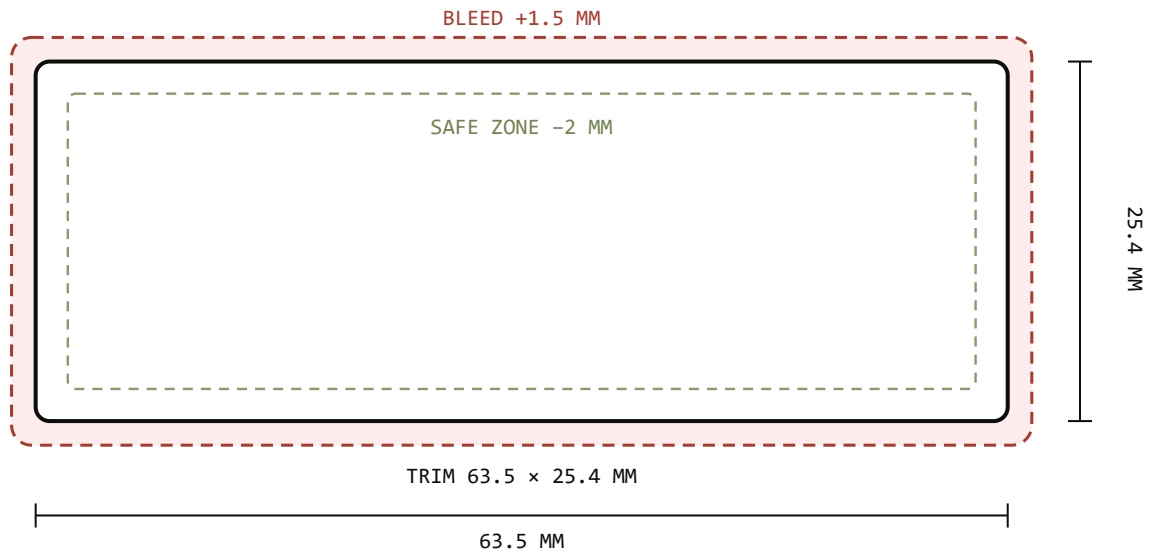


Dieline spec

- Trim: **63.5 × 25.4 mm**
- Bleed: **+1.5 mm** all sides (full-bleed art)
- Safe zone: **-2 mm** inside trim - all text and QR inside
- Corner radius: **1.5–2 mm**

File setup

- Artwork at 100% scale, 300 dpi minimum for raster
- Vector logo and type; fonts outlined on export
- Export PDF/X-1a or PDF/X-4, CMYK, no RGB objects
- Dieline on its own spot-color layer, not printing

The centered stack, dimensioned.

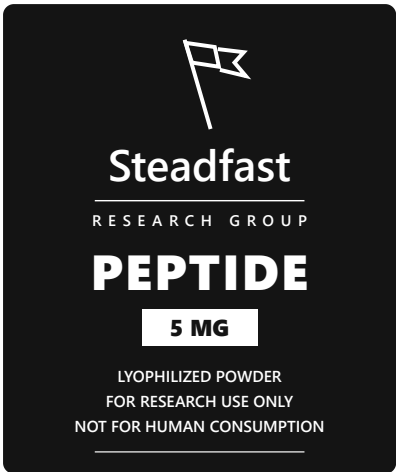

#	ELEMENT	SPEC AT 3 ML SCALE
1	Flag mark	White outline, ≈8 mm tall, centered
2	Steadfast	Manrope SemiBold, mixed case, ≈4.2 mm cap height
3	Rule	0.4 pt white, ≈80% of panel width
4	RESEARCH GROUP	Caps, +0.35em tracking, ≈1.6 mm
5	Product name	Manrope ExtraBold caps, ≈4.8 mm (PEPTIDE = placeholder per SKU)
6	Strength box	White rectangle; text in the label color; ≥5.5 pt
7	Compliance	3 centered lines, caps, ≥5 pt, never removed
8	Rule	Closes the stack

COLORWAY	LABEL	CAP	ELEMENTS
01 • Black	#141414 matte	Black flip-off	All white; strength text black
02 • Sage	≈#869274 matte	Sage flip-off	All white; strength text sage

- Everything sits on one centered vertical axis on the front panel (≈40% of wrap width, opposite the seam).
- Lot, expiry, and the COA QR print on the back of the wrap - they never crowd the front panel.
- Per SKU, only two things change: the product name and the strength box. Nothing else moves.

Hierarchy rule: at arm's length the read order is product name → strength → brand → compliance. If anything competes with the name, shrink it.

Built for a freezer, printed like a flagship.

Color at print

COLOR	HEX	CMYK (APPROX.)	USE ON LABEL
 Carbon	#0D0D0D	60 / 40 / 40 / 100	Background floods only
 White	#FFFFFF	0 / 0 / 0 / 0	Type on Carbon
 Ink	#23271F	0 / 0 / 0 / 100 (text)	Type on light labels - plain 100K, no rich black
 Olive	#8A9A6B	47 / 26 / 70 / 5 • ≈PMS 5777	Verified data, accents, cap match
 Bone	#F4F1E8	3 / 3 / 8 / 0	Light-system label stock tint

- Small text is always plain 100K (or pure white reversed). Rich black on 5 pt type mis-registers and blurs.
- Check contrast: olive on Carbon passes for accents only - never for compliance text.
- For the flag printed small (<8 mm), use the solid one-color version; the outline fills in.

Material stack

LAYER	SPEC	WHY
Facestock	White BOPP (polypropylene), 50–60 µm	Waterproof, tear-proof; paper wrinkles with condensation
Finish	Matte laminate (preferred) or matte varnish	Premium feel, no glare under lab light, protects print
Adhesive	Permanent acrylic, freezer-grade (apply ≥10°C, service to –40°C)	Vials live at –20°C and sweat on thaw
Print	Digital (short run) or flexo (≥5k); thermal-transfer overprint for lot/expiry	Variable data stays crisp and smudge-proof

The unboxing is part of the standard.

PIECE	SIZE (W × D × H)	CONSTRUCTION	FINISH
Single-vial carton	30 × 30 × 60 mm	Reverse-tuck end, 350 gsm SBS board, optional vial retainer flap	Matte Carbon outside, olive interior flood, white flag + wordmark
10-vial kit box	120 × 65 × 60 mm	Rigid two-piece or telescope lid; die-cut foam or molded pulp insert, 10 × Ø23 mm cavities	Matte Carbon, foil-stamped or spot-UV flag, mono batch panel on base
Insert card	70 × 100 mm	400 gsm uncoated, single card	Reconstitution + storage guidance, QR to COA and education library
Shipper	Fits kit + 2 gel packs	Plain corrugate + insulated liner for summer lanes	Unbranded exterior (discretion), flag stamp inside lid

Packaging rules

- Compliance line appears on every printed piece: carton base panel, kit box base, insert footer.
- Lot and expiry are printed (not stickered) on the carton end flap at ≥6 pt.
- Exterior shipper stays unbranded; the brand moment happens on opening.
- One accent maximum per piece - olive flood *or* foil flag, never both.

Sequence on open: Carbon box → olive interior → vial with the flag facing up. Check this on the packing line; it is specified, not accidental.

Before anything goes to print.

- ☐ Physical vial from the actual glass lot measured; label size confirmed against it
- ☐ Trim, bleed (+1.5 mm), and safe zone (–2 mm) set; dieline on non-printing layer
- ☐ Peptide name is the dominant element; hierarchy matches Section 03
- ☐ Compliance line present, mono, ≥5 pt: FOR RESEARCH USE ONLY · NOT FOR HUMAN CONSUMPTION
- ☐ Lot + expiry zone reserved for thermal overprint; QR ≥10×10 mm, tested from 30 cm
- ☐ QR resolves to the correct batch COA URL (https, no redirect chain)
- ☐ Small text 100K plain black or pure white; no rich-black type
- ☐ Olive matched to cap sample under daylight, not screen
- ☐ Fonts outlined; PDF/X export; 300 dpi effective resolution on raster
- ☐ Freezer-grade adhesive + matte laminate confirmed on the print order
- ☐ One printed proof approved on the actual vial before the full run

Regulatory note: research-use-only peptide labeling still carries legal exposure. Have counsel review final label text (claims, disclaimers, state requirements) before the first production run. This guide covers design specification, not legal compliance.